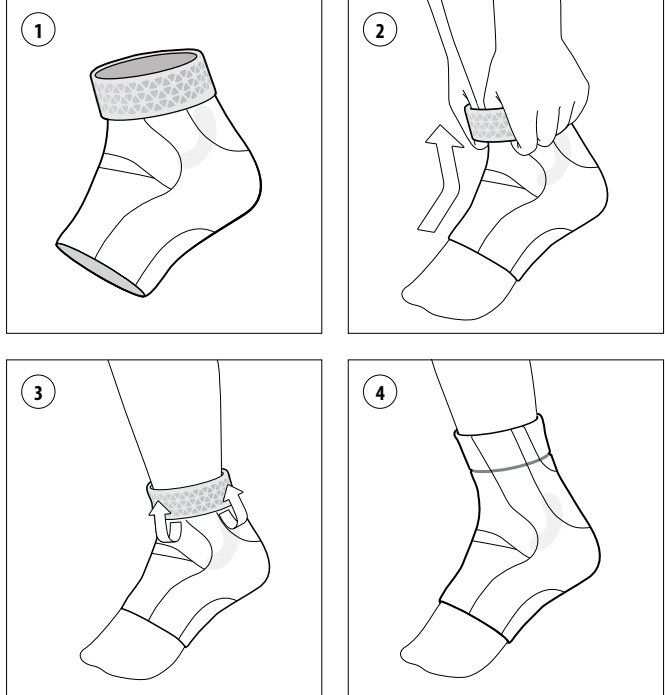
ENGLISH

SYMBOLS

MD	Medical Device
-----------	----------------

LATEX	Caution: This device contains natural rubber latex which may cause allergic reactions.
--------------	--

INTENDED USE

The device is intended for increased proprioception and compression of the ankle

Indications for use

For Ankle conditions that may benefit from compression, such as;

- Ankle injury rehabilitation
- Mild sprains and strains
- Mild instability

Contraindications

Must not be used by individuals for whom compression is contraindicated.

Warnings and Cautions:

- The device should not be used for extended periods of time and removed during long periods of rest, e.g., sleep.
- Patients with already compromised blood circulation in their extremities, like those with diabetes or peripheral vascular disease, need to be extra careful and are advised to consult with a healthcare professional before using the device. If you experience pain, swelling, sensation changes, or if your extremity shows signs of insufficient blood flow (turns blue, white, or cold) while using this device, discontinue the use immediately and consult with a healthcare professional.
- The device should be used directly on intact skin and never over an open wound.

GENERAL SAFETY INSTRUCTIONS

Read these instructions carefully before use. Keep them for future reference.

Any serious incident in relation to the device must be reported to the manufacturer and relevant authorities.

The patient should stop using the device and contact a healthcare professional:

- If there is a change or loss in device functionality, or if the device shows signs of damage or wear hindering its normal functions.
- If any pain, skin irritation, or unusual reaction occurs with the use of the device.

The device is for single patient – multiple use.

USAGE

Size Selection

Measure the ankle circumference.

The sizing chart is shown on the package, If between sizes always select the smaller size for optimal fit.

Device Application

- Fold the top part of the device outward so that the FoldLock edge becomes visible (**Fig. 1**)
- Grasp the device with two hands on the side and pull it over the foot until it reaches the heel. (**Fig. 2**)
- Unfold the top part of the device so that the FoldLock edge comes in contact with the skin (**Fig. 3**).

Device removal

Fold the top part of the device outward so that the FoldLock edge becomes visible. Grip the device on the side and push it off in a downward direction.

Cleaning and care

Wash regularly to retain the compression of the textiles.

- Machine-wash with a mild detergent at cold temperature (30 °C)
- Air dry.

Note: Secure the hook and loop fasteners before washing to avoid damaging the textiles.

Note: Do not tumble dry, iron, bleach, or wash with fabric softener.

Note: Avoid contact with salt water or chlorinated water. In case of contact, rinse with fresh water and air dry.

DISPOSAL

The device and packaging must be disposed of in accordance with respective local or national environmental regulations.

LIABILITY

Össur does not assume liability for the following:

- Device not maintained as instructed by the instructions for use.
- Device assembled with components from other manufacturers.
- Device used outside of recommended use condition, application, or environment.

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