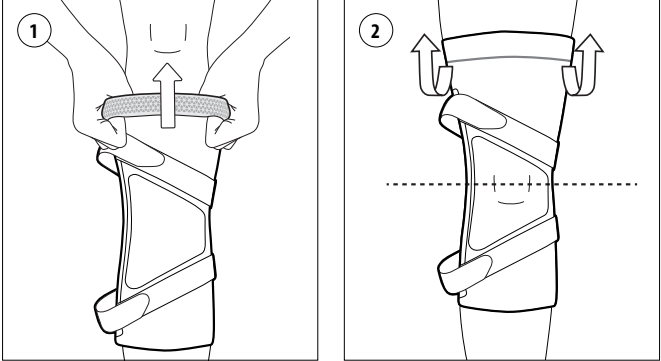
- For easy donning, fold the top part of the support outward so that the FoldLock edge becomes visible. Grasp the support with two hands on the side and pull the support over the knee with the hinge on the affected side of the knee (**Fig. 1**).
- Pull evenly over the knee joint until the knee cap is free from pressure and centered into the patella graphics. Then fold the edge of the support back into the original position (**Fig. 2**).
- Bend the knee to an 80° angle and tighten the DFS straps until pain relief is sufficient (**Fig. 3**). With the knee in full extension, verify correct position and strap tension. Readjust, if needed (**Fig. 4**).

Device removal

Fold the FoldLock edge of the support outward so that the edge becomes visible. Undo the DFS straps. Attach the hook and loop straps appropriately, to avoid damaging the knit. Grip the support at the level of the calf and pull it off in a downward direction.

Cleaning and care

Wash regularly to retain the compression of the textiles.

- Machine-wash with a mild detergent at cold temperature (30°C)
- Air dry

Note: Secure the hook and loop fasteners before washing to avoid damaging the textiles.

Note: Do not tumble dry, iron, bleach, or wash with fabric softener.

Note: If the device is used in salt water or chlorinated water, rinse with fresh water and air dry.

DISPOSAL

The device and packaging must be disposed of in accordance with respective local or national environmental regulations.

LIABILITY

Össur does not assume liability for the following:

- Device not maintained as instructed by the instructions for use.
- Device assembled with components from other manufacturers.
- Device used outside of recommended use condition, application, or environment.

ENGLISH

SYMBOLS

	Medical Device
	Caution: This device contains natural rubber latex which may cause allergic reactions.

INTENDED USE

The device is intended for unicompartmental unloading of the knee

Indications for use

For osteoarthritis and other knee conditions that may benefit from unloading and/or compression, such as:

- Medial or lateral knee pain
- Early or mild uni-compartmental osteoarthritis of the knee
- Symptomatic, uni-compartmental meniscus degeneration

Contraindications

Must not be used by individuals for whom compression is contraindicated.

Warnings and Cautions:

- The device should not be used for extended periods of time and removed during long periods of rest, e.g., sleep.
- Patients with already compromised blood circulation in their extremities, like those with diabetes or peripheral vascular disease, need to be extra careful and are advised to consult with a healthcare professional before using the device. If you experience pain, swelling, sensation changes, or if your extremity shows signs of insufficient blood flow (turns blue, white, or cold) while using this device, discontinue the use immediately and consult with a healthcare professional.

GENERAL SAFETY INSTRUCTIONS

Read these instructions carefully before use. Keep them for future reference.

Any serious incident in relation to the device must be reported to the manufacturer and relevant authorities.

Never wear the device directly over an open wound. Inspect the skin under the device's contact area before and after each use.

The patient should stop using the device and contact a healthcare professional:

- If there is a change or loss in device functionality, or if the device shows signs of damage or wear hindering its normal functions.
- If any pain, skin irritation, or unusual reaction occurs with the use of the device.

The device is for single patient – multiple use.

USAGE

Size Selection

Measure the calf and thigh circumference 15 cm (6”) above and below mid-patella. The sizing chart is on the packaging, if between sizes always select the smaller size for optimal fit.

Device Application

- First loosen the two ends of the DFS straps without pulling them out of the slits.

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