
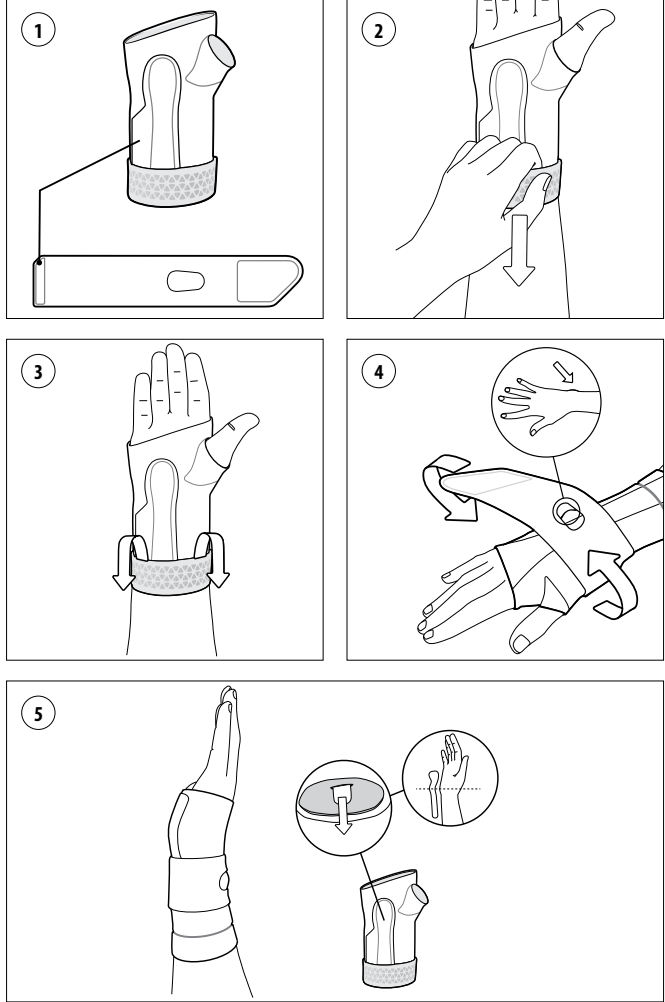
ENGLISH

SYMBOLS

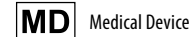INTENDED USE

The device is intended to provide support and compression of the wrist

Indications for use

For wrist conditions that may benefit from compression, such as:

- Rheumatic Arthritis (RA) and other chronic pain conditions.
- Wrist sprain and strains.
- Overuse and pain conditions.
- Wrist tendinitis.

Contraindications

Must not be used by individuals for whom compression is contraindicated.

Warnings and Cautions:

- The device should not be used for extended periods of time and removed during long periods of rest, e.g., sleep.
- Patients with already compromised blood circulation in their extremities, like those with diabetes or peripheral vascular disease, need to be extra careful and are advised to consult with a healthcare professional before using the device. If you experience pain, swelling, sensation changes, or if your extremity shows signs of insufficient blood flow (turns blue, white, or cold) while using this device, discontinue the use immediately and consult with a healthcare professional.
- The device should be used directly on intact skin and never over an open wound.

GENERAL SAFETY INSTRUCTIONS

Read these instructions carefully before use. Keep them for future reference.

Any serious incident in relation to the device must be reported to the manufacturer and relevant authorities.

Warning: If there is a change or loss in device functionality, or if the device shows signs of damage or wear hindering its normal functions, the patient should stop using the device and contact a healthcare professional. The device is for single patient – multiple use.

USAGE

Size Selection

Measure wrist circumference at the level of the radial styloid with the hand outstretched. Wrist in

0° of flexion. The sizing chart is shown on the package, if you are between sizes always select the smaller size for optimal fit.

Device Application

- Fold the top edge of the device outward so that the FoldLock edge becomes visible (**Fig 1**).
- Grasp the device and pull over the hand, so that the wrist is enclosed (**Fig 2**).
- Unfold the top part of the device so that the FoldLock edge is in contact with the skin (**Fig 3**).
- Fasten the RadialStrap onto the device and wrap it around the wrist. The bony prominence on the outside of the wrist (Ulnar Styloid Process) should be enclosed by the oval cut-out (**Fig 4**).
- The aluminum stay on the palm of the hand can be removed from the device in order to bend to suit the contours of your wrist. To do this, pull the stay out of the pocket on the inside of the device (**Fig 5**).

Device removal

Undo the RadialStrap and fold the top edge of the device outward so that the FoldLock edge becomes visible. Grip the device and pull it off the wrist.

Cleaning and care

- Machine-wash with a mild detergent at cold temperature (30°C)
- Air dry

Wash regularly to retain the compression of the textiles.

Note: Do not tumble dry, iron, bleach, or wash with fabric softener.

Note: Secure the hook and loop fasteners before washing to avoid damaging the textiles.

Note: Avoid contact with salt water or chlorinated water. In case of contact, rinse with fresh water and air dry.

DISPOSAL

The device and packaging must be disposed of in accordance with respective local or national environmental regulations.


LIABILITY

Össur does not assume liability for the following:

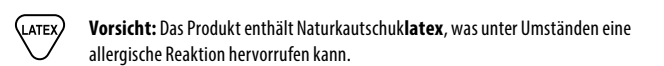- Device not maintained as instructed by the instructions for use.
- Device assembled with components from other manufacturers.
- Device used outside of recommended use condition, application, or environment.

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